

POWERFUL SUCTION

Well dried, the vac-clean uv 2[®] is used with a powerful vacuum suction capacity of over 60 kPa. For this purpose, the attachment head optimized for the type of cleaning is attached to the basic head. At the start, the vacuum suction power with sonic vibrations ensures a powerful cleaning effect. The attachment heads can be easily guided around openings, crevices, bushes or in a gentle circular motion. The base head collects all dirt without reducing the suction power. The suction power can be varied in 3 modes at the touch of a button. The blue LED light (450nm) illuminates the area to be cleaned and also has an antimicrobial hygiene function. The base head and the attachment heads can be easily changed and cleaned in an ultrasonic bath (e.g. u-sonic 3).

The new vac-clean uv 2[®] is a medical device and meets the high requirements for medical devices (Class I) of the Medical Device Regulation (EU) 2017/745.

vac-clean uv 2[®]

Rechargeable mini vacuum cleaner for hearing aids and hearables.




Your specialist

- ✓ Powerful vacuum pump
- ✓ 7 attachment heads for different applications
- ✓ Medical device according to MDR (EU)

VERSATILE

The new rechargeable mini vacuum cleaner vac-clean uv 2[®] is suitable for all hearing systems (BTE, ITE, RIC, cochlear) all charging stations, all OTC hearing aids, all True Wireless Stereo (TWS), a variety of electronic products such as cell phones, remote controls and much more. These products suffer from adhering dirt and moisture. Failures are the result. The remedy is continuous cleaning and care.

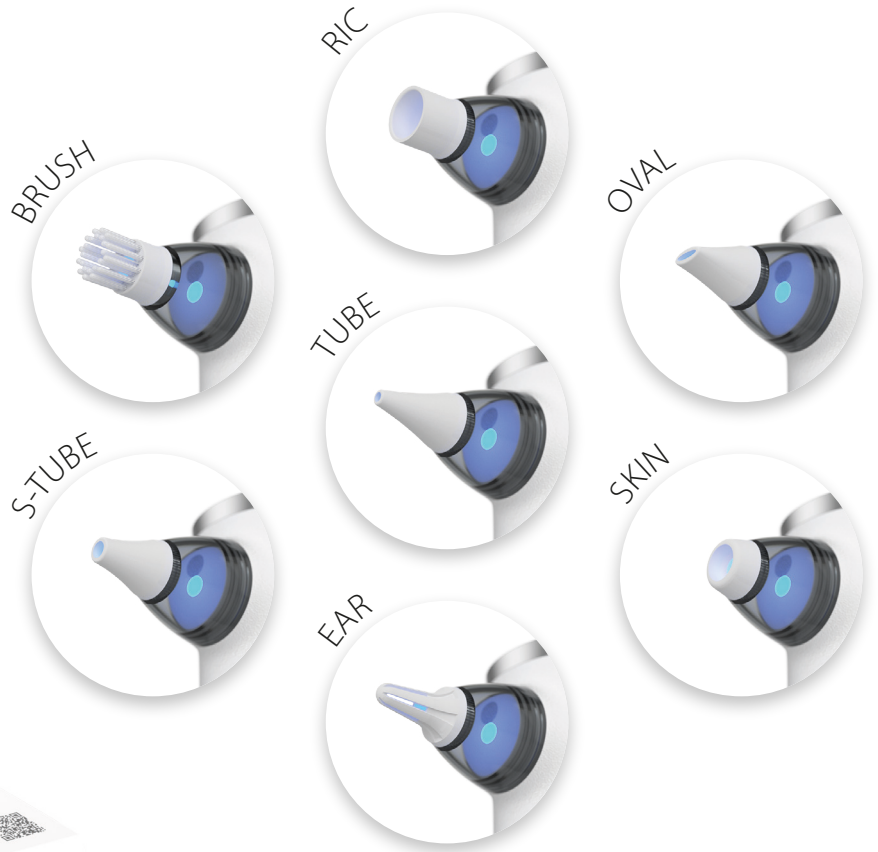

Effective suction by vacuum pump



Sonic-Brush for perfect cleansing



Blue LED light for illumination



Easy handling by only one tactile button



Rechargeable by LiPo battery



Mobility by USB-C standard



3 intensity levels for flexible application



All features controlled by microprocessor



Highest product quality with 2 year warranty



Energy-efficient, gentle and absolutely safe